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INDEX

PRESENTATION OF MANUAL	24
IDENTIFICATION OF EQUIPMENT	25
Product Identification	25
Principles and fundamentals applied to the operation of the product	26
Description of Equipment	26
Indication of the equipment	26
MODULES, ACCESSORIES, OPTIONS AND MATERIALS OF CONSUMPTION.....	27
TECHNICAL SPECIFICATIONS	28
- Technical characteristics	28
- Symbolologies of packaging	29
- Symbolologies of product	29
- Contents of the accessible and inaccessible markings	30
EQUIPMENT OPERATION	32
- Conditions for installation and use	32
- Handling container of the dosage	33
- Control and operation panel	34
- Selection memory	35
PRECAUTIONS, RESTRICTIONS AND WARNINGS	37
- Conditions of transport and storage	37
- Sensitiveness to environmental conditions foreseeable in normal situations of use ..	37
- Precautions and warnings "during the installation" of equipment	38
- Recommendations for preserving the equipment.....	38
- Precautions and warnings "during the use" of equipment	38
- Precautions and warnings "after" the use of equipment	38
- Precautions and warnings during the "cleaning and disinfection" of equipment	39
- Precautions to be adopted against foreseeable or uncommon risks, related to the deactivation and abandoning of equipment.....	39
CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION.....	39
- Additional procedures for reuse	39
- Cleaning and disinfection.....	39
- Preventive maintenance	41
- Corrective maintenance	41
UNFORESEEN EVENTS – SOLUTION OF PROBLEMS	41
WARRANTY OF EQUIPMENT	42
FINAL CONSIDERATIONS.....	42

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"Declared exempt from registration by the Ministry of Health"

ATTENTION

For greater safety:

Read and understand all the instructions contained in these Instructions for Use before installing or operating this Equipment.

Note: These Instructions for Use must be read by all the operators of this Equipment.

IDENTIFICATION OF THE EQUIPMENT

Dear Customer

Congratulations. You have made a good choice when you decided to buy a GNATUS QUALITY product comparable to the best products available in the World. This manual is a general presentation of your product and it will give you important details to help you to solve possible problems.

Please, read it and keep this with you.

Identification

Technical name: Alginate Mixer

Brand: GNATUS

Nome Comercial: Algimixer

IDENTIFICATION OF THE EQUIPMENT

Description of the Equipment

Automatic mixer for handling alginate to prepare dental impressions. Compact and simple to use, it is composed by 4 containers, 2 spatulas, stirrer, liquid container and a plastic basket with magnet (for fixation).

The Algimixer provides time reduction in the mixture preparation due to its constant speed. With memorization for 3 different times, the device can be used on up to 16 different speeds (digital display), with beep at the time of completion of the procedure, ensuring a uniform mixture, free of air bubbles. The equipment also possesses an automatic lock that offers a safe procedure in the process of alginate Manipulation.

Principles and fundamentals applied to the operation of the product

The Automatic Alginate mixer is composed by an electric engine with stable movement, it has a digital display for viewing time and memorization.

It is also consisted of 4 containers, 2 spatulas, stirrer, 1 liquid container and a plastic basket with magnet (for fixation).

Indication of Equipment

This equipment must be used and handled by a trained person (properly regulated professional, according to local legislation of country) observing the instructions contained in this manual.

It is the user's obligation to use the equipment only in perfect condition and protect himself / herself, patients and others from possible dangers.

MODULES, ACCESSORIES, OPTIONS AND MATERIALS OF CONSUMPTION


The content of this page is for information purposes, and the device may be presented different from that illustrated. Therefore, when purchasing the product check the technical compatibility between the device, coupling and accessories.

The drawing illustrates the equipment with all the optional items, your equipment will consist only of the items chosen for its purchase option.

The use of any part, accessory or material not specified or referred to in these instructions for use is the user's responsibility.

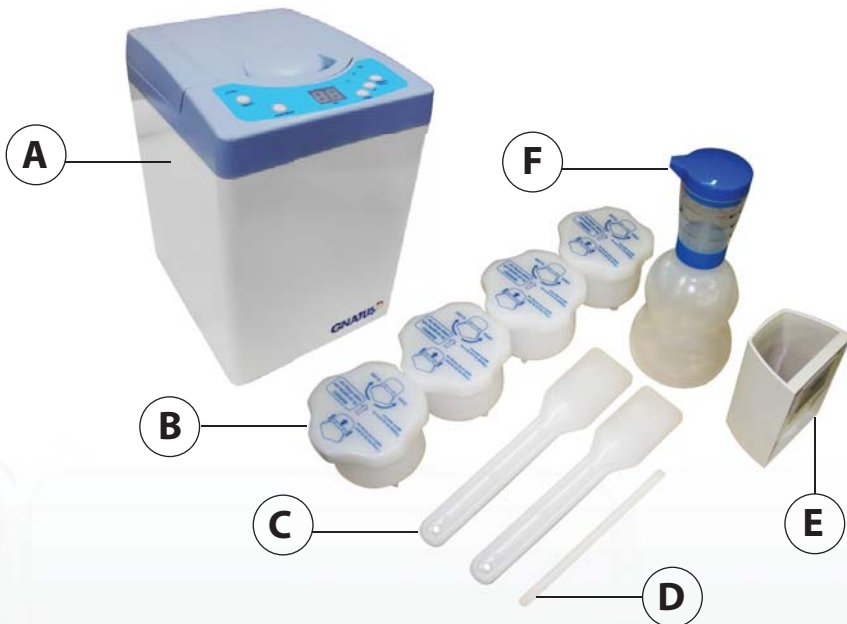- A - Alginate Mixer
- B - Container
- C - Spatula
- D - Stirrer
- E - Plastic basket
- F - Bottle for liquids

TECHNICAL SPECIFICATIONS

Technical characteristics of equipment and accessories

General

Model:
Alginate Mixer

Power

Power Supply:
127V ou 240V
Frequency:
50/60Hz
No. of stages:
Monophasic
Power:
350W
Rated Current:
127V - 4,4A(50Hz) and 3A(60Hz)
220V - 1,4A(50Hz) and 1,5A(60Hz)

Other Specifications

Dimensions:
29,5mm x 245mm x 20,5mm
Net weight:
16,4 Kg
Gross weight: :
18 kg
R.P.M máxima:
3600

Resources:

Memory Select:
Selects 3 programming options (editable), pre-configured with 8, 10 and 12 sec
Timer:
1-16 seconds - Adjustable by second

TECHNICAL SPECIFICATIONS

The materials used to produce the equipment are Biocompatible.

The parts of this equipment are not subject to degradation by wear, corrosion, fatigue or aging of materials during their lifetime.

Parts lists and circuit diagrams

Gnatus states that the provision of circuit diagrams, parts list or any other information that will provide technical assistance by the user, can be requested, if previously agreed between the Company and Gnatus.

Package symbols

Maximum stacking: It determines the maximum quantity of boxes which can be stacked during transportation and storage "as per packaging".

Packing to be transported and / or stored avoiding humidity, rains and wet floor.


Packing to be transported and / or stored with the harrows up.

The packing must be stored and transported away from direct sun light exposure.

Packing to be transported and / or stored with care (should not suffer drop and neither receive impact).

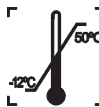Temperature limit for the packing to be stored or transported.

Product symbols

Careful: It indicates an important instruction for the operation of the product. Not following it can cause dangerous malfunctioning.

Important: It indicates an instruction of safety for operation of the product. Not following it, can lead to serious danger to the patient.

Note: It indicates useful information for operation of the product.

Warning - see the manual

TECHNICAL SPECIFICATIONS

Product symbols

Source Voltage

Type B

Display digital

Time setting button

Start/Stop

General key

Material recyclable electronics

Grounding (at several points of the equipment) indicates the condition of being grounded.

On/Off

Button time setting

Button temperature setting

Memory select

Fuse

Content of accessible and non-accessible demarcations

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EQUIPMENT OPERATION

Conditions for installation and use:

1 - Remove the screws

- When removing the Algimixer from the box, remember that you must remove the two protection screws, located below. **Fig. A.**
- Insert the protection screws behind the Algimixer, in order to keep them stored. **Fig. B.**

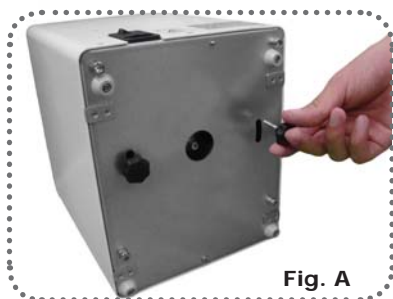
Fig. A

Fig. B

Important: These two protection screws are used whenever the Algimixer device is transported, with the function of protecting the internal components.

2 - Electrical power

- Insert the power cable in the power input plug. **Fig. C.**
- Press the main switch on the back of the product. **Fig. D.**

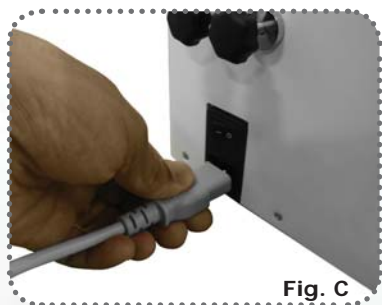
Fig. C

Fig. D

Check the voltage of the equipment at the time of making the electric installation.

EQUIPMENT OPERATION

Handling of the dosage container

- Unscrew the bottom cover as **fig. A**.
- When inserting water, close the cover.
- The dispenser contains measure points with 3 dosage levels, **fig. B**.
- When pressing the bottle, you control the amount of water that will be inserted into the dispenser according to the Alginate dosage.
- To obtain higher viscosity, use 3 Alginate doses in only 2 levels of water, **fig. C**.

Fig. A

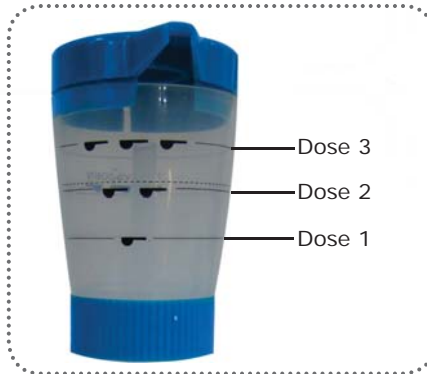
Fig. B


Fig. C

English

Prepare the mixture according to the instructions of the alginate supplier.

- To open / close the mixing container, just follow the procedure of figures below:

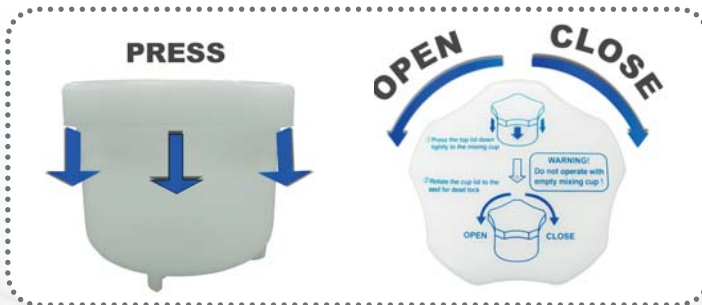Attention: Do not operate the equipment with empty mixing container!

EQUIPMENT OPERATION

Control and operations panel

- Fit the container inside the Algimixer.

Attention: Never remove the protection band **(A)**.

Make sure the container is correctly closed so it will not open during the mixing process. Do not start the procedure with the cover opened.

A - ON / OFF - Turn on and off the panel

B - START / STOP - Start and stop.

C - SEC - Time Viewer.

D - DOWN - Reduce the time.

E - UP - Increase the time.

F - MEMORY SELECT - Selection of memorized time.

G - Memory time indicator LED related to dosage.

EQUIPMENT OPERATION

Memory Select

1- Press the "memory select" button to select the time required for the mixture. The "memory select" option is preset by the manufacturer with 3 different options

Memory	Mixing Time (sec.)
1	8
2	10
3	12

2- If you want to change the memory of the mixing time established by the manufacturer, you can make the adjustment by pressing the "memory select" to select 1 of 3 times to be edited, then pressing "UP" or "DOWN" to enter the new time to be memorized. The maximum time is 16 seconds. You can adjust it second by second. To complete click (hold down) on the "MEMORY SELECT" button until hearing an audible signal (beep).

3- The mixing time should be adjusted according to the listed instruction pattern;

- For each dose, a longer time is needed in Algimixer.
- The time setting for the power network frequency of 50Hz (water temperature between 10 ° C ~ 20 ° C) should proceed with the following measurement as below:

Dose	Mixing Time (sec.)
1	7 ~ 10
2	8 ~ 11
3	9 ~ 12

- The time setting for the power network frequency of 60Hz (water temperature between 10 ° C ~ 20 ° C) should proceed with the following measurement as below:

Dose	Mixing Time (sec.)
1	6 ~ 9
2	7 ~ 10
3	8 ~ 11

EQUIPMENT OPERATION

Memory Select

If there is any difference between the water temperature and the room temperature established by the alginate supplier, the setting time will be different.

Cleaning

After use, the debris inside the container must be removed by performing a cleaning.

The small hole in the center of the cover shall be carefully cleaned, otherwise, the cover

will not be opened or closed easily.

If there is dirt in the Algmixer, use a clean cloth, or use a neutral cleaning agent.

Do not use volatile solvent oil for cleaning, it may wash-out the equipment surface.

Replacing the fuses

- Before exchanging fuses, disconnect the plug from the power network.
- With the aid of a screwdriver, loosen the cover (**Fig. A**) of the fuse holder, and then replace the fuse by the spare fuse (**Fig.B**).

Note: The spare fuse accompanies the equipment, after making the first exchange it is advisable to always have a backup fuse.

Fig. A

Fig. B

PRECAUTIONS, RESTRICTIONS AND WARNINGS

Transportation and storage

This equipment must be transported and stored observing the following directions:

- Avoid falls and impacts;
- Keep it dry, do not expose it to rain, water drops or wet floor;
- Keep it away from water and direct sunlight, and in its original wrapping;
- Don't move it over irregular surfaces, protect it from rain and observe the maximum stack quantity specified in the packaging;
- Transportation and storage temperature range: -12°C to 50°C;
- Transportation and storage relative humidity range: 0°C to 90°C;
- Use the screws to protect the region below the equipment, not to damage the inside (Fig.A - Pag.10).

The Equipment maintains its condition of safety and efficacy, provided that it is maintained (stored) as mentioned in this instruction of use. Thus, the equipment will not lose or alter its physical and dimensional features.

Sensitivity to environmental conditions in normal situations of use

- The equipment has been planned not to be sensitive to interference such as magnetic fields, external electrical factors, electrostatic discharge, pressure or variance of pressure, provided that the equipment is installed, maintained, clean, preserved, transported and operated as per this instruction for use.

Precautions and warnings “during the installation” of equipment

- Make sure that the voltage of the equipment is in accordance with the local network (127 or 220V).
- Position the equipment in a place where it will not get wet.
- Install the equipment in a place where it will not be damaged by pressure, temperature, humidity, direct sunlight, dust, salts, or sulfur compounds.
- This equipment was not planned for use in an environment where vapors, anesthetic mixtures inflammable with air, or oxygen and nitrous oxide can be detected.
- Check the voltage of the equipment at the moment of executing the electrical installation.
- The equipment must be grounded correctly.
- Before the first use and/or after long interruptions from work such as vacations, clean and disinfect the equipment; eliminate air and water deposited in the internal hoses.

Grounding:

The grounding is required, only in this way the user is in complete safety. It is the responsibility only to the user by not grounding, implying the loss of warranty.
Use three-prong outlet. Never connect the neutral grounding.

PRECAUTIONS, RESTRICTIONS AND WARNINGS

Recommendations for the dental equipment maintenance

Your Gnatus equipment has been designed and developed according to the standards of modern technology. Similarly to other kinds of equipment, it requires special care, which is many times neglected due to several reasons and circumstances.

Therefore, here are some important reminders for your daily routine. Try to follow these simple rules, which will save you a lot of time and will avoid unnecessary expenses once they start making part of your working procedure.

Precautions and warnings “during the use” of equipment

- The equipment should only be operated by duly enabled and trained technicians (Dental Surgeons, Capacitated Professionals)
- If any maintenance should be required, only use services of the Gnatus Authorized Technical Assistance.
- The equipment has been manufactured to handle both continuous and intermittent operation; so follow the cycles described in these Instructions for Use.
- Although this equipment has been planned in accordance with the standards of electromagnetic compatibility, it can, in very extreme conditions, cause interference with other equipment. Do not use this equipment together with other devices very sensitive to interference or with devices which create high electromagnetic disturbance.
- Do not expose the plastic parts to contact with chemical substances, use in the routines of dental treatment, such as: acids, mercury, acrylic liquids, amalgams, etc.
- Never block the ventilation openings of the device. Keep openings free of lint, hair and other materials.
- Do not use the product without filling the tank with water. The dry running can damage the unit.
- Do not fill the tank above the maximum to avoid overflow.

Gnatus shall not be responsible for:

- Use of the equipment differing from that for which it is intended.
- Damages caused to the equipment, the professional and/or the patient by the incorrect installation and erroneous procedures of maintenance, differing from those described in these Instructions for use which come with the equipment or by the incorrect operation of it.

Precautions and warnings “after” the use of equipment

- Turn off the main switch of the dental set when it is not in use for an extended period of time.
- Always maintain the equipment clean for the next operation.
- Do not modify any part of the equipment. Do not disconnect the cable or other connections without need.
- After using the equipment, clean and disinfect all the parts which may be in contact with the patient.

PRECAUTIONS, RESTRICTIONS AND WARNINGS

Precautions and warnings during the “cleaning and disinfection” of equipment

- Before cleaning the equipment, turn off the main switch.
- Avoid spilling water, or other liquids inside the equipment, which could cause short circuits.
- Do not use microabrasive material or steel wool when cleaning, or employ organic solvents or detergents which contain solvents such as ether, stain remover, gasoline etc.
- When cleaning of silver or copper, where the effect of oxidation caused a blackout, you must add special solutions to the water to remove rust.

Precautions in case of alteration in the functioning of equipment

- If the equipment has any abnormality, check if the problem is related to any item listed in the topic of unforeseen events (failures, causes and solutions). If it is not possible to resolve the problem, turn off the equipment, remove the power supply cable from the socket and contact your representative (Gnatus).

Precautions to be adopted against foreseeable or uncommon risks, related to the deactivation and abandoning of equipment

In order to avoid environmental contamination or undue use of the Equipment after it has become useless, it should be discarded in the suitable place (as per the local legislation of the country).

- Pay attention to the local legislation of the country for the conditions of installation and disposal of residue.

CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION

Additional procedures for reuse

The equipment can be reused in undetermined, i.e. unlimited, quantities, only needing to be cleaned and disinfected.

Cleaning and disinfection

Important: In order to execute cleaning or any type of maintenance, ensure that the equipment is disconnected from the electrical network.

The cleaning procedure below should be made to start work and after each patient.

Always turn off the main switch before performing the daily maintenance procedures

CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION

After the exchange of fluid in the stainless steel bowl, remove the waste deposited in the bottom by using a soft cloth or sponge. Dry and clean by using a cloth moistened with alcohol 70°.

To clean the body, including the cover, use a cloth moistened with neutral soap and water or the product recommended below:

To clean the equipment, we recommend the use of **“BactSpray”** (Reg n° MS: 3.2079.0041.001-5) or any other similar product:

Active component:
Benzalkonium chloride (tri-quaternary ammonium)
Solution 50%..... 0.329%
Chemical composition: Butyl Glycol, Decyl polyglucose, Sodium Benzoate, Sodium Nitrate, Essence, Deodorized Propane / Butane, demineralized Water.

For more information concerning cleaning procedures, see manufacturer’s instructions.

ATTENTION:

- This product can also be used for cleaning and disinfection of the water unit basin.
- In order to eliminate safety risks or damage to equipment, we recommend avoiding liquids from getting inside the equipment when cleaning.
- The application of other solvent-based or sodium hypochlorite chemical products to clean is not recommended, since they can damage the equipment.

 **NOTE:** The registration at the Ministry of Health of the “BactSpray” is executed separately from the product described in this manual, as the “BactSpray” is not manufactured by Gnatus.

 **NOTE:** Use gloves and other systems of protection, during the disinfection.

CORRECTIVE AND PREVENTIVE MAINTENANCE AND PRESERVATION

Preventive Maintenance

In order to protect your equipment, seek Gnatus technical assistance for periodic revisions of preventive maintenance.

Corrective Maintenance

If the equipment has any abnormality, check if the problem is related to any of the items listed in the item Unforeseen Events (situation, cause and solution). If it is not possible to solve the problem, turn off the equipment, and request Gnatus technical assistance.

UNFORESEEN EVENTS – SOLUTION OF PROBLEMS

⚠ Upon coming across any problem in operation, follow the instructions below to check and repair the problem, and/or get in touch with your representative.

Problem	Probable cause	Solution
-Inoperative apparatus	<ul style="list-style-type: none"> - Disconnected plug from the outlet. - Blown fuse. - Lack of energy in the network. -The cover is not totally closed 	<ul style="list-style-type: none"> -Connect the plug. -Check the fuse and exchange it, if necessary. -Wait for network normalization. - Close the mixer cover
- Display off (no image).	<ul style="list-style-type: none"> - Faulty board - Faulty driving engine or gear. 	<ul style="list-style-type: none"> -Request the presence of a Gnatus authorized technician. -Request the presence of a Gnatus authorized technician.
- Slow startup of apparatus	<ul style="list-style-type: none"> - Faulty driving engine or gear. - Faulty startup capacitor 	<ul style="list-style-type: none"> -Request the presence of a Gnatus authorized technician. -Request the presence of a Gnatus authorized technician.
-Abnormal noise or vibration during operationo	<ul style="list-style-type: none"> - Excess of Alginate and/or water -Empty container -Damage to the gear -Alginate dust spill because of loose cover 	<ul style="list-style-type: none"> -Maximum number of mixtures 100 grams. -Never use the mixer with the empty container -Request the presence of a Gnatus authorized technician. -Check if the container cover is closed tight.

WARRANTY OF EQUIPMENT

This equipment is covered by the warranty terms counting from the date of installation, as specified below; provided that the defect has occurred in normal conditions of use and that the equipment has not remained stored for more than 06 months counting from the issue date of the sales document until the date of the actual installation.

- WARRANTY TERMS: TWELVE MONTHS;

Note: It is not covered by this warranty the following item: Fuses.

- LOSS OF THE WARRANTY:

A) Attempt to repair using an inadequate tool or by unauthorized technicians;

B) Installation of the equipment by an unauthorized technician;

C) Damage arising from inappropriate storage or signs of infringement;

D) Incorrect use of the equipment;

E) Use of a cleaning product not indicated by the factory;

F) Falls or blows which the equipment may undergo or lack of observation of an compliance

with the guidelines of the Owner's Manual, which was delivered with the present document, together with the equipment. Repair or replacement of parts during the warranty period shall not extend the validity term of their warranty.

- This warranty does not exempt the customer from paying the service charge for the visit and the travel expenses of the technician, except when the customer sends the equipment to execute the maintenance inside the establishment of the technical assistance.

"Consumer Defense Code - art. 50, unique paragraph".

- The Warranty Certificate comes with the product and must be filled in upon the date of installation by the Gnatius Authorized Technician.

- Queries and information: GNATUS Help Desk (+55) 16 2102-5000.

FINAL CONSIDERATIONS

The most important aspect related to equipment care is that concerning spare parts.

To guarantee the life span of your equipment, use only original Gnatius spare parts. They are sure to follow the technical specifications and standards required by Gnatius.

We must also point out to you our chain of authorized dealers. Only dealers that make part of this chain will be able to keep your equipment constantly new for they count on technical assistants who have been trained and on specific tools for the correct maintenance of your equipment.

Doubts and information: GNATUS Call center (55-16) 2102-5000.

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